

**Important Note:** The anatomical position of the Nalu Neuromodulation System implanted in the patient must be known prior to performing the MRI examination. The supervising physician (e.g., the radiologist) must review the information provided by the implanting physician, the patient's Medical Device Identification Card, and/or obtain an x-ray to determine the anatomic location of the implanted Nalu Neuromodulation System in the patient's body.


Additionally, the patient must undergo a proper MRI screening procedure according to the MRI facility's policies and procedures.

**NEUROLOGY/REFERRING PHYSICIAN:** Please complete the PATIENT ELIGIBILITY FORM and include information pertaining to the location of the implanted Nalu Neuromodulation System. Provide the PATIENT ELIGIBILITY FORM to your patient to take to the MRI Center and make a copy to include in the patient's medical file.

**MRI/RADIOLOGY:** Patient safety will be ensured by using the information on this form and filling out the PRE-MRI CHECKLIST to verify the use of the acceptable conditions, with special attention to the location of the Nalu Neuromodulation System and in consideration of the anatomic area of interest for the MRI examination. Additionally, specific conditions must be followed to ensure patient safety relative to the prevention of excessive heating of the Nalu Neuromodulation System.

**Important Note:** In addition to utilizing the PATIENT ELIGIBILITY FORM AND PRE-MRI CHECKLIST, carefully review the latest Nalu Neuromodulation System Instructions for Use (IFU) available online at [nalumed.com](http://nalumed.com) and follow all relevant information before proceeding with the MRI examination.

### *MRI Safety Information*

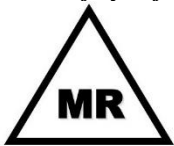

### **MR Conditional**

**WARNING:** An MRI examination performed outside the guidelines may result in the MRI-related electromagnetic fields adversely interacting with this implanted device, potentially injuring the patient and damaging the device.

**NALU MEDICAL TECHNICAL AND CLINICAL SUPPORT IS AVAILABLE MONDAY THROUGH FRIDAY, 8 AM – 5 PM PACIFIC STANDARD TIME PLEASE CALL US AT (800) 618-3402 or visit [www.nalumed.com](http://www.nalumed.com)**

**PATIENT ELIGIBILITY FORM:** To be filled out by the referring physician and/or appropriate Nalu personnel prior to MRI.

Date: \_\_\_/\_\_\_/\_\_\_

### CONTACT INFORMATION

Patient Name (last name/first name): \_\_\_\_\_


Referring Physician (last name/first name): \_\_\_\_\_

Physician's Phone #: \_\_\_\_\_

### NALU NEUROMODULATION SYSTEM

Model Number(s): \_\_\_\_\_

Nalu Neuromodulation System Position: \_\_\_\_\_

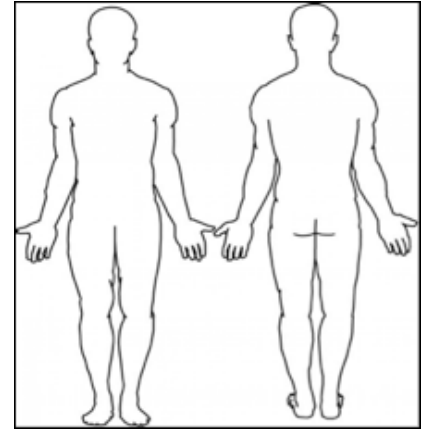

*Sketch the position of the Nalu Neuromodulation System*

**VALID IMPEDANCE CHECK:** Can be completed by the patient using the Nalu R/C application or appropriate Nalu personnel prior to MRI.

Date: \_\_\_\_\_

Pass / Fail: \_\_\_\_\_

**WARNING:** This impedance measurement is considered valid for seven (7) days for the purposes of an MRI scan as long as the patient has not experienced any trauma to the implanted system or significant change in related pain or paresthesia.

**All of the components of the Nalu Neurostimulation System are MR Unsafe except for the following:**

- *Nalu Implantable Pulse Generators (Dual 8-Contact Integrated Model 11002-040, Single 8-Contact Integrated Model 11001-040)*
- *Nalu Implantable Pulse Generators (Single 4-Contact Integrated 25 cm Tined Model 11005-025, Single 4-Contact Integrated 40 cm Tined Model 11005-040)*
- *Nalu Implantable Pulse Generators (Dual 8-Contact Ported Model 11004-002 connected to one or two 40 cm Nalu Lead Model 12001-040 or 12007-040, Single 8-Contact Ported Model 11003-002 connected to a 40 cm Nalu Lead Model 12001-040 or 12007-040)*
- *Nalu Implantable Pulse Generator (Single 4-Contact Ported Model 11006-002 connected to a 25 cm Nalu Tined Lead Model 12005-025 or a 40 cm Nalu Tined Lead Model 12005-040, Dual 4-Contact Ported Model 11007-002 connected to one or two 25 cm Nalu Tined Lead Model 12005-025 or one or two 40 cm Nalu Tined Lead Model 12005-040 or one 40 cm Nalu Tined Lead Model 12005-040 and one 25 cm Nalu Tined Lead Model 12005-025)*
- *Nalu Anchor (Model 13001)*




**PATIENT ELIGIBILITY FORM AND PRE-MRI CHECKLIST**

- *All implantable components from convenience kits Model 71001, Model 71002, Model 71003, Model 71004, Model 71005, Model 71006, Model 71007, MODEL 71010, Model 71019, Model 71020, Model 71021, Model 71022, Model 72001, Model 72006, Model 72007, and Model 72011*

**MEASURING THE IMPEDANCE**

To measure the impedance of electrodes:

1. Tap on **Check Impedance** on the menu.
2. Tap on the **Measure Impedance** on the Check Impedance screen.
3. Tap on the **View Result** button. **Electrode Impedance** screen appears. It displays the most recent Impedance values with date and time.

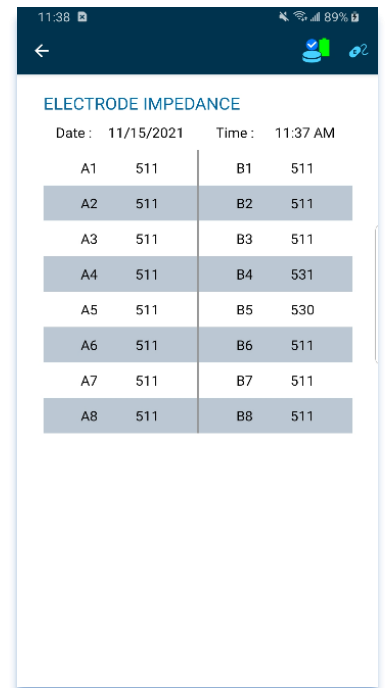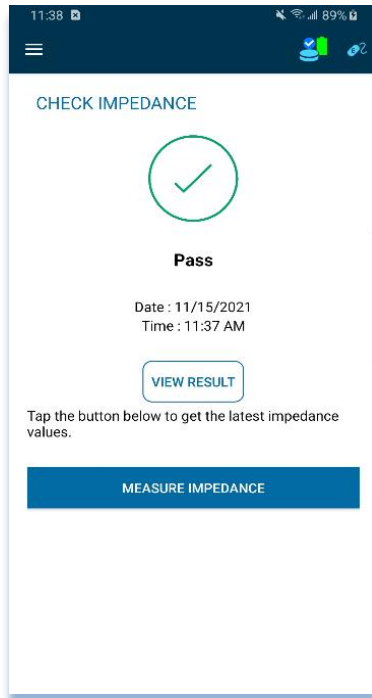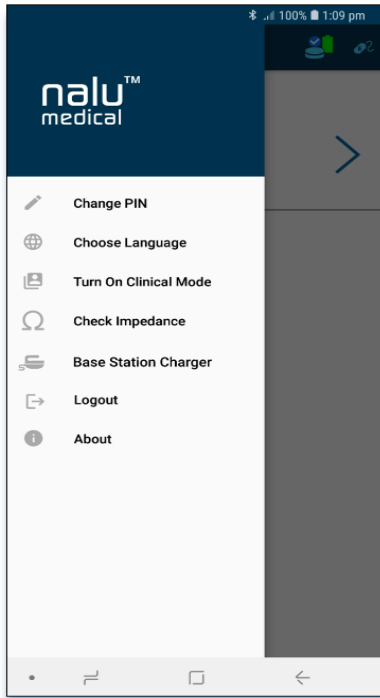

**NOTE: You can view the most recent impedance values for up to 48 hours.**

Impedance results can be **Pass**, **Fail** and **Inconclusive** based on the impedance values of the configured electrodes.

Impedance Result	Condition
Pass	If Impedances of the configured electrodes are in-range (0 to 9500 $\Omega$ ).
	If the number of leads is unspecified by Remote Control and in-range impedances (0-9500 $\Omega$ ) are detected for electrodes.
Fail	If Impedance of any configured electrodes is above 9500 $\Omega$ .
Inconclusive	If the number of leads is unspecified by Remote Control and in-range impedances (0-9500 $\Omega$ ) are not detected for electrodes.

**NOTE:** You cannot measure impedance if the IPG/EIC is disconnected before or during an impedance measurement. A notification saying **“Please connect IPG/EIC to measure impedance”** appears.

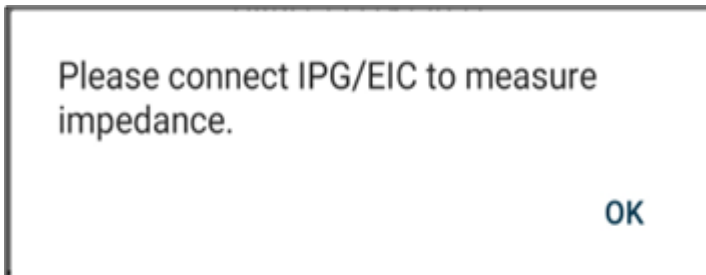


Figure 11 Notification message

**CHECKLIST:** To be filled out by MRI/Radiology personnel prior to MRI.

The MRI examination will be performed on the patient with the implanted Nalu Neuromodulation System because there is a valid indication as determined by the supervising physician and all guidelines must be carefully followed to ensure patient safety.

**Preparation of the Patient Prior to the MRI Examination**

- Inform the patient of the risks associated with undergoing an MRI examination: an MRI exam performed outside recommended guidelines may result in the electromagnetic fields used with MRI technology interacting adversely with an implanted Nalu Neurostimulation System, potentially injuring the patient and/or damaging the device.
- A trained healthcare professional with the proper knowledge of MRI technology such as an MRI safety-trained radiologist, MRI technologist, MRI nurse, or MRI physicist must ensure that the MRI examination will be conducted according to the information presented in this document.
- Document the patient’s programming parameters.
- Perform an impedance check within 7 days of the MRI examination. Do not perform an MRI if the impedance is greater than 9500 Ω.
- Do not sedate or anesthetize the patient so that the patient can inform the MRI system operator of any unusual sensations or problems associated with the MRI examination.
- Instruct the patient to immediately inform the MRI system operator if any discomfort, stimulation, shocking, or heating is experienced during MRI.
- All external components of the Nalu Neuromodulation System are MR Unsafe and are contraindicated for the MR system room. Therefore, the Therapy Disc, Wearable Garment, and Adhesive Clip must be removed before the patient is allowed into the MR system room.

Therapy Disc:	<input type="checkbox"/> Removed	<input type="checkbox"/> N/A Note: _____
Adhesive Clip:	<input type="checkbox"/> Removed	<input type="checkbox"/> N/A Note: _____
Wearable Garment:	<input type="checkbox"/> Removed	<input type="checkbox"/> N/A Note: _____

**Additional Preparation of Patients with more than one Nalu Neurostimulation System**

- If the patient has more than one implanted Nalu Neurostimulation System, there must be a minimum distance of 50 mm between any part of the systems (IPG and lead(s)). Do not perform an MRI examination on a patient if there is less than 50 mm between the implanted systems.
- If the Nalu Neurostimulation Systems are implanted in different body regions, it can be assumed that there is at least 50 mm between the systems. The body regions, as referred to above, are defined as the head, neck, torso, upper extremities, and lower extremities.
- If the Nalu Neurostimulation Systems are implanted in the same body region, the distance between the systems must be confirmed before the MRI examination. The distance should be confirmed through imaging, such as fluoroscopy or X-Ray.

**Additional conditions for all MRI Examinations**

- Do not perform an MRI if the patient has lead extension or leads from a different manufacturer attached to the Nalu Implantable Pulse Generator. The risk of performing an MRI examination under those circumstances has not been evaluated and, thus, may cause harm to the patient and/or the components.
- Do not perform MRI on a patient undergoing the trial phase of the Nalu Implantable Pulse Generator (i.e., the patient has a percutaneously implanted lead and an external Trial Therapy Disc (Model 34002)).
- Do not perform MRI on a patient that has any other active medical implants.

**Spinal Cord Stimulation (SCS) or Peripheral Nerve Stimulation (PNS)**

**Head and Extremities**

**Each condition indicated below must be carefully followed and verified to ensure patient safety:**

- Static magnetic field of 1.5 Tesla (T)/64 MHz or 3.0 T/128 MHz, only.
- Maximum spatial field gradient of 2,000 gauss/cm (20 T/m).
- Maximum MRI system reported, whole body averaged specific absorption rate (SAR) of 2.0 W/kg for 15 minutes of scanning (i.e., per pulse sequence) in the Normal Operating Mode.
- For head/brain MRI examinations, only use a transmit/receive RF head coil. **No parts of the implanted Nalu Neurostimulation System may be within the transmit/receive RF head coil.**
- For extremity MRI examinations, only use a transmit/receive RF coil that includes a knee, foot/ankle, or wrist transmit/receive RF coil. **No part of the implanted Nalu Neurostimulation System is allowed within the transmit/receive RF coil used.**
- A valid impedance measurement has been performed within 7 days of the MRI examination.

Under the scan conditions defined above, the Nalu Neurostimulation System is expected to produce a maximum temperature rise of 2.5°C after 15 minutes of continuous scanning (i.e., per pulse sequence) when using one of the above types of transmit/receive RF coils at a maximum, whole body averaged SAR of 2 W/kg. In non-clinical testing, the image artifact caused by the Nalu Neurostimulation System extends approximately 10 mm from this implant when imaged using a gradient echo pulse sequence and a 3 T MRI system.

**Spinal Cord Stimulation (SCS) or Peripheral Nerve Stimulation (PNS)**

**Full Body**

**Each condition indicated below must be carefully followed and verified to ensure patient safety:**

- Use the manual MA-000105 entitled, “Magnetic Resonance Imaging (MRI) Safety Information,” which pertains to performing an MRI examination safely on a patient. The manual identifies the appropriate MRI labeling conditions for the patient’s implanted IPG and lead, nerve target, IPG location, and the orientation of the lead’s stimulating contacts.
- Static magnetic field of 1.5 Tesla (T)/64 MHz only.
- Maximum spatial field gradient of 2,000 gauss/cm (20 T/m).
- Up to 30-minute exposure to whole-body continuous scanning with a waiting period of 30 minutes before the next imaging session. The SAR value shall adhere to the RF conditions listed in the MRI Safety Information table on the patient’s corresponding label.
- A valid impedance measurement has been performed within 7 days of the MRI examination .

Under the scan conditions defined above, the Nalu Neurostimulation System is expected to produce a maximum temperature rise of 6°C after 30 minutes of continuous scanning (i.e., per pulse sequence) when using full body and the RF conditions specified on the patient’s corresponding label. In non-clinical testing, the image artifact caused by the Nalu Neurostimulation System extends approximately 10 mm from the implant when imaged using a gradient echo pulse sequence and a 3 T MRI system.

**Considerations during the MRI Examination**

- Similar to other MRI examinations, carefully monitor the patient throughout the MRI procedure both visually and audibly. Immediately discontinue the MRI examination if the patient reports any problems or unusual sensations.

**Considerations after the MRI Examination**

- After the patient leaves the MRI system room, turn the Therapy Disc on and verify connection to the Implantable Pulse Generator.
- Perform an impedance check.